


Given by Markus Seeberger over Teams 17.10.2025

Setting Up QUARC Software

Software Sources

QUARC Real-Time RapidControl: <https://www.quanser.com/products/quarc-real-time-control-software/>

Our software license covers up to version 4.2.65535.65535 (Quarc 2022 SP1).

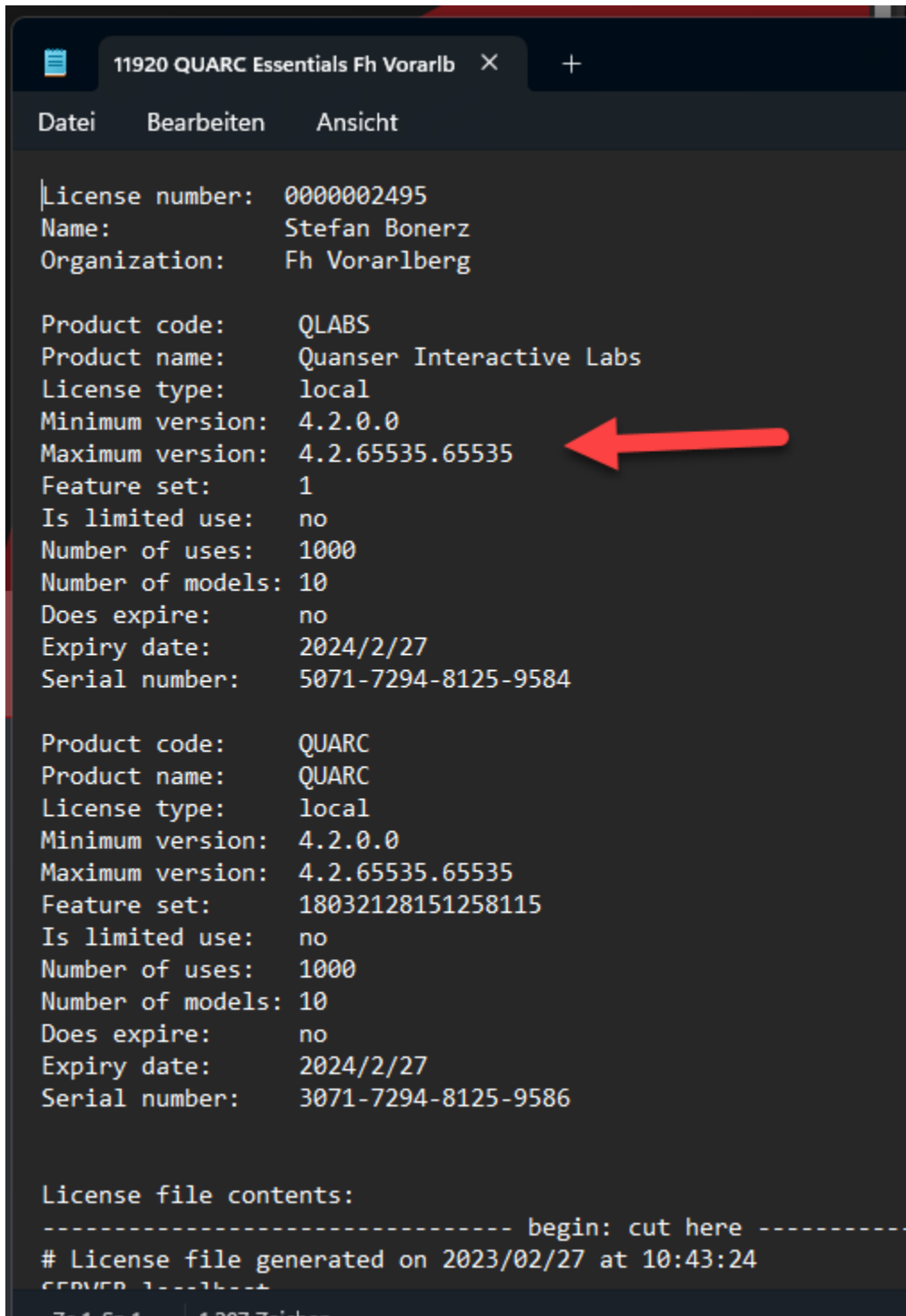

```
11920 QUARC Essentials Fh Vorarlb X +
Datei Bearbeiten Ansicht

License number: 0000002495
Name: Stefan Bonerz
Organization: Fh Vorarlberg

Product code: QLABS
Product name: Quanser Interactive Labs
License type: local
Minimum version: 4.2.0.0
Maximum version: 4.2.65535.65535
Feature set: 1
Is limited use: no
Number of uses: 1000
Number of models: 10
Does expire: no
Expiry date: 2024/2/27
Serial number: 5071-7294-8125-9584

Product code: QUARC
Product name: QUARC
License type: local
Minimum version: 4.2.0.0
Maximum version: 4.2.65535.65535
Feature set: 18032128151258115
Is limited use: no
Number of uses: 1000
Number of models: 10
Does expire: no
Expiry date: 2024/2/27
Serial number: 3071-7294-8125-9586

License file contents:
----- begin: cut here -----
# License file generated on 2023/02/27 at 10:43:24
C:\Program Files\Quanser\Quanser Interactive Labs\
Zeile 50: 1.207 Zeichen
```


NI Compact RIO (Maximum version 18.0 or 19.6):

<https://www.ni.com/de/support/downloads/drivers/download.ni-compactrio.html#>

Visual Studio Community: <https://visualstudio.microsoft.com/de/vs/community/>

Installation (in this exact order)

1. Install Matlab with Simulink. Matlab must be version 2022b at most. For newer versions, a newer Quarc license is required. Simply install all components. License code is in the file FIK.txt and the network license file is network.lic.

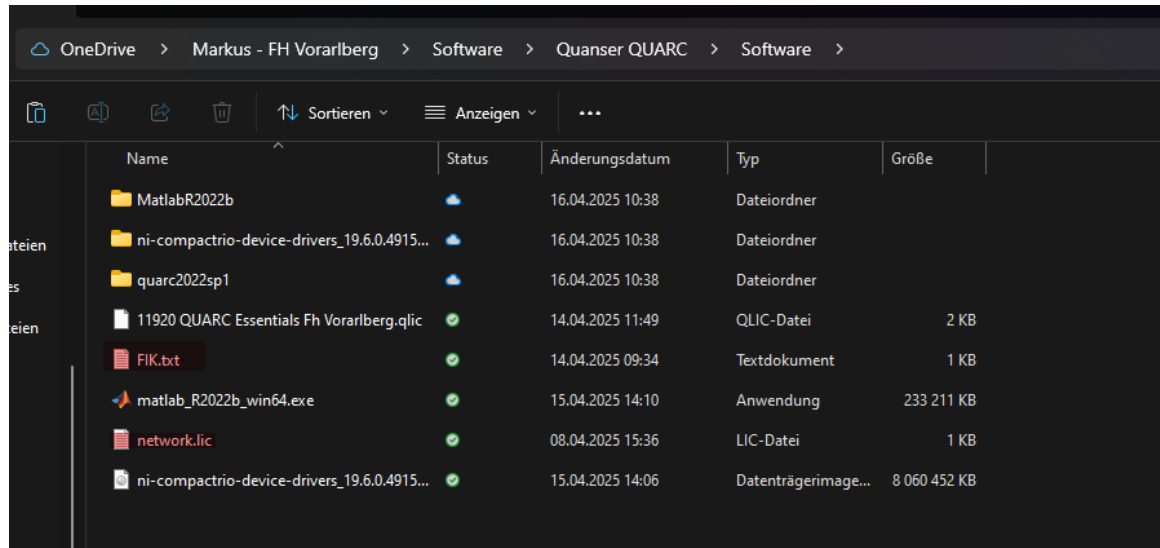

- 2.
2. Install Visual Studio Community for C++ Desktop (Version 17 was used).
3. Run the command 'mex -setup' in the command line for Matlab integration.



```
Command Prompt
Microsoft Windows [Version 10.0.26100.3775]
(c) Microsoft Corporation. All rights reserved.

C:\Users\Admin>mex -setup
MEX configured to use 'Microsoft Visual C++ 2022 (C)' for C language compilation.

To choose a different language, execute one from the following:
mex -setup C++
mex -setup FORTRAN


C:\Users\Admin>
```

4. Install NI Compact RIO (Maximum version 18.0 or 19.6). Install completely.
5. Install QUARC Real-Time RapidControl Software (Quarc 2022 SP1). If the installation order was not followed, please perform the steps under 'Quarc_library not found' to register the library.

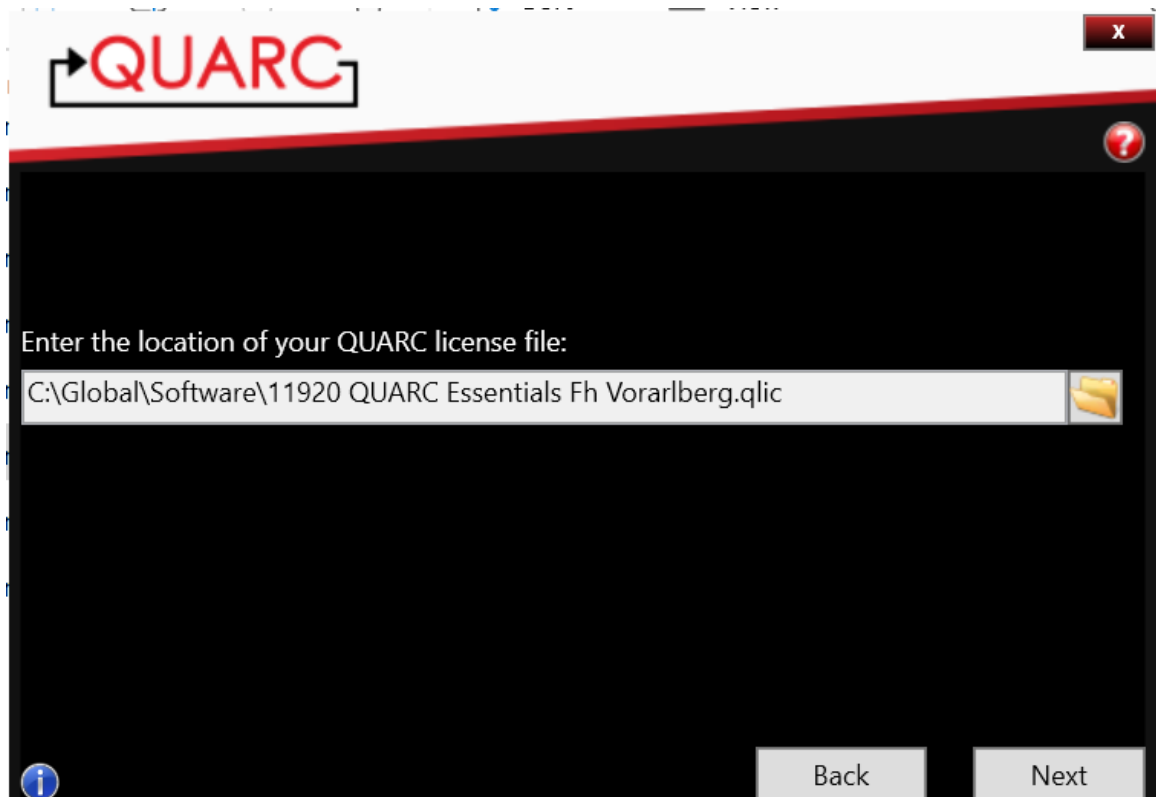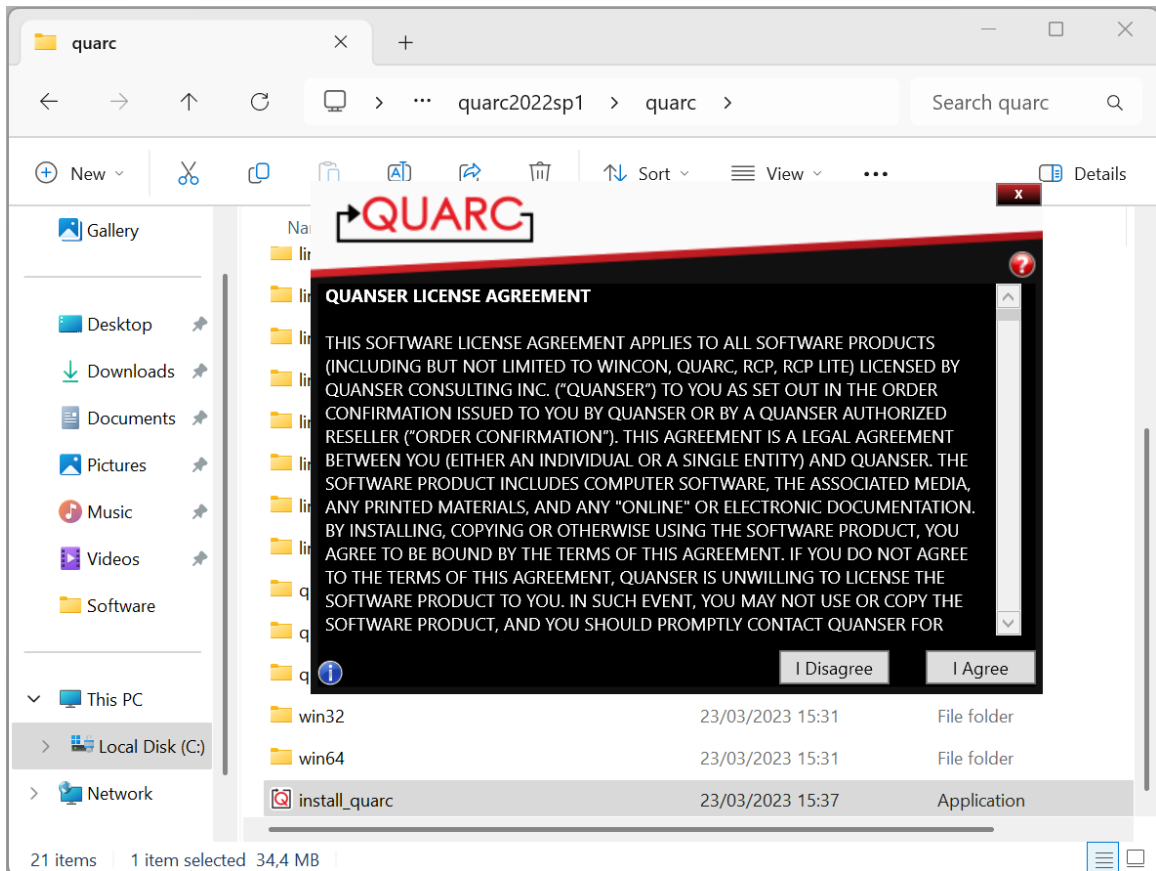


Given by Markus Seeberger over Teams 17.10.2025

If the above order is followed, all dependencies should be fulfilled.

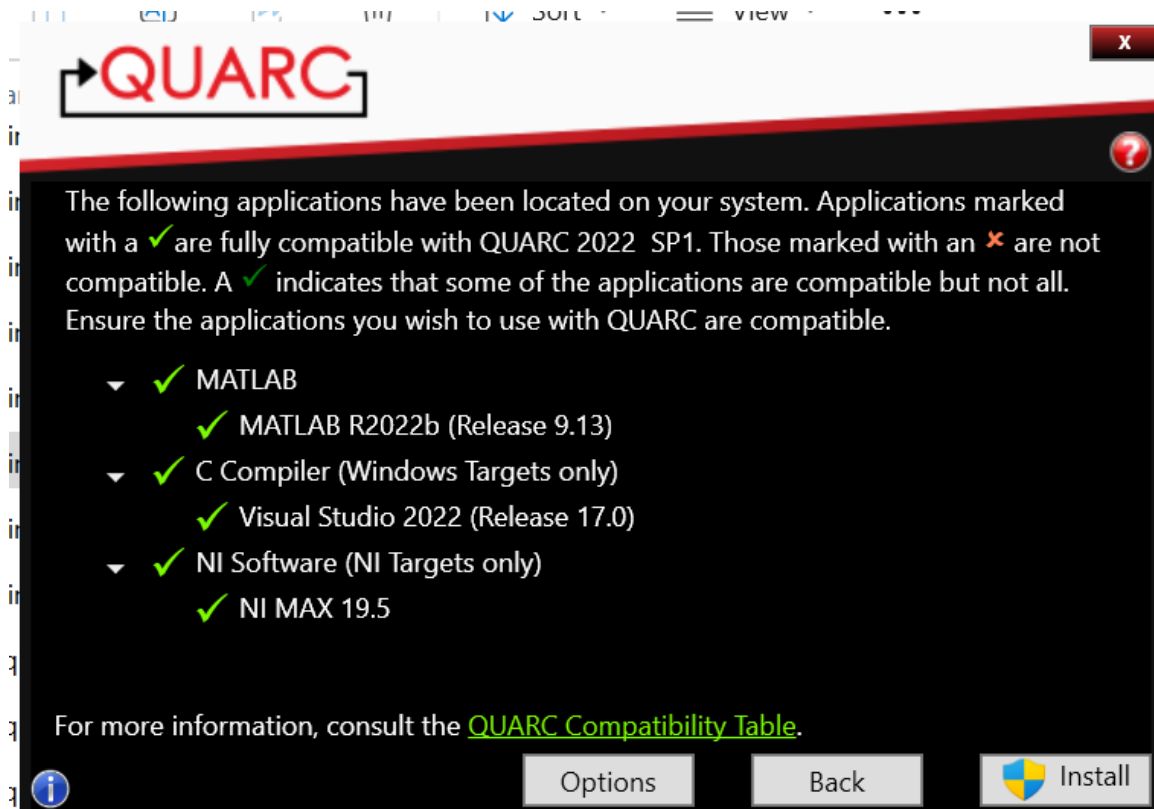


After successful installation, the included sample project in Simulink should look like this:

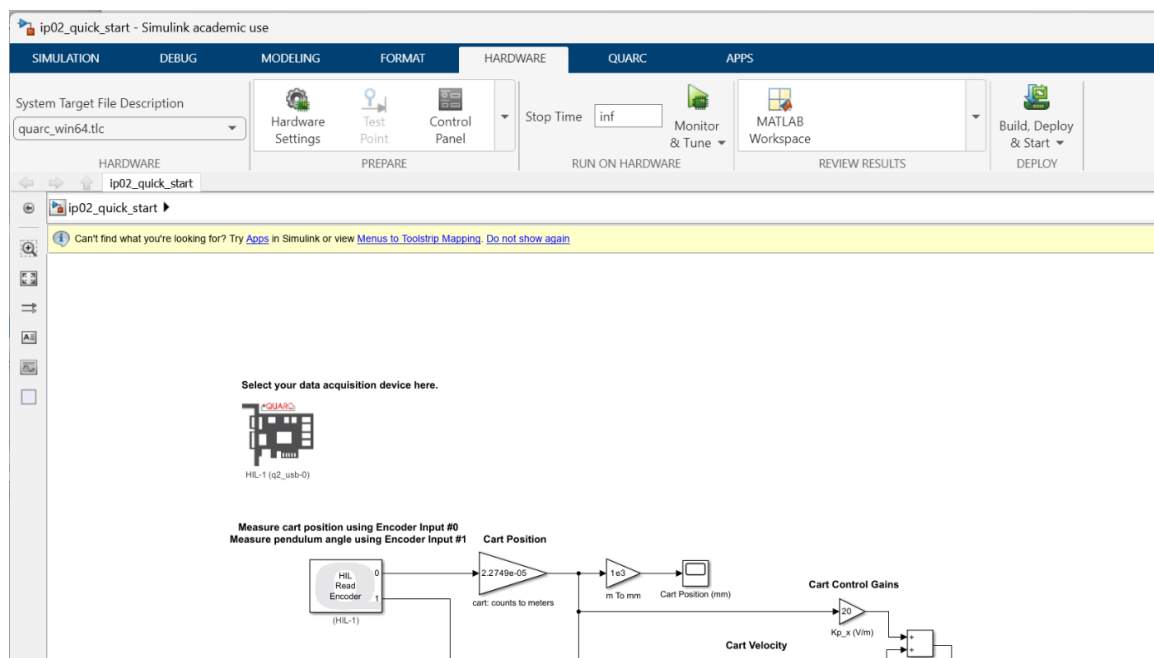

Troubleshooting

'Quarc_library' not found

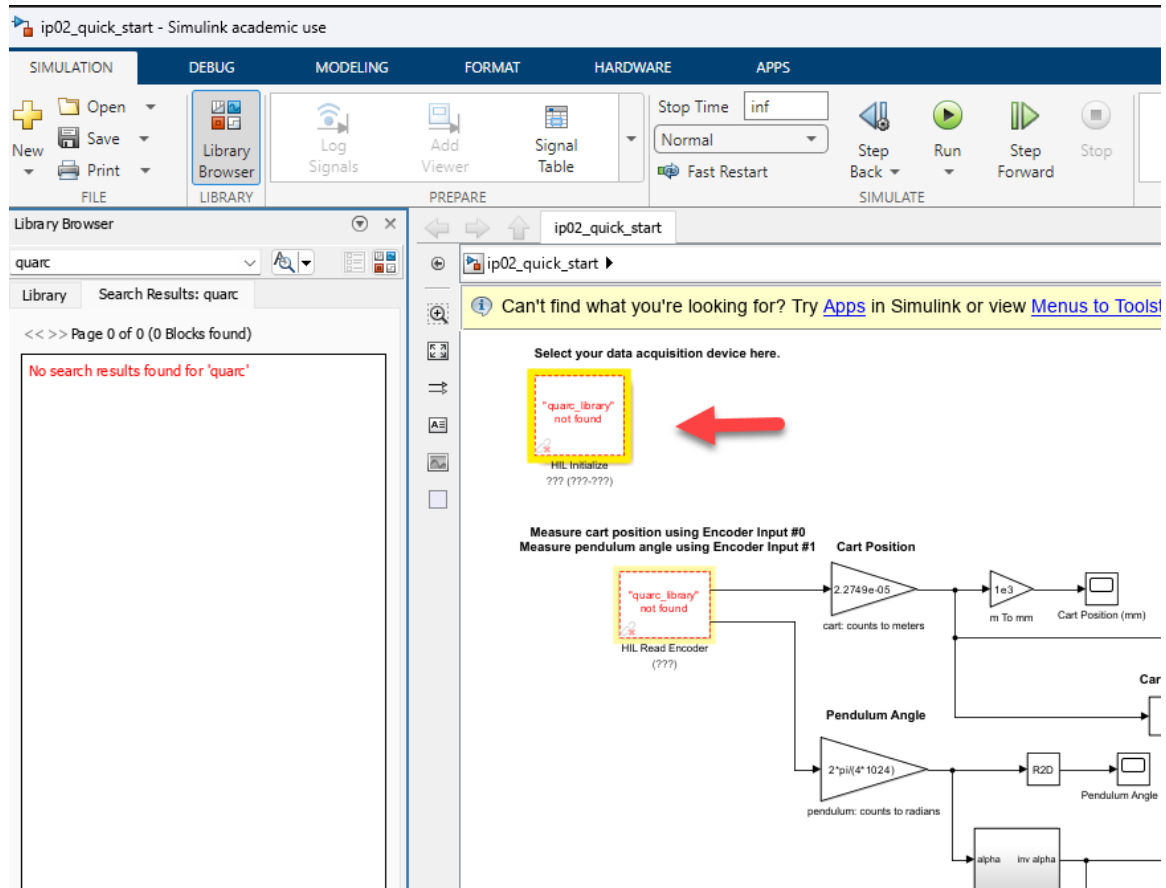

Solution according to: https://www.quanser.com/faq/im-getting-a-quarc_library-not-found-error-in-my-simulink-blocks/

The QUARC Real-Time Control software may not be installed or QUARC may have been installed before MATLAB. If QUARC was not installed, then close MATLAB and install the QUARC Real-Time Control software as instructed in the Quick Installation Guide. If QUARC was installed and MATLAB was installed after, then run the following commands:

```
cd(fullfile(getenv('QUARC_DIR'), 'quarc'));
```

```
quarc_setup;
```

This function adds the QUARC Targets library to the MATLAB path and updates the MATLAB documentation. It also performs other setup tasks to make QUARC fully functional in MATLAB and Simulink.

Notes by Thomas Lønne Stiansen

Toolchain not recognized / Installation Repair / License

Upon major updates of Windows, the Quanser programs may run into trouble. This makes it unable to recognize the toolchain selected in Simulink, which is solved by running the repair installation: **C:\Global\Software\quarc2022sp1\quarc\install_quarc.exe**

This repair has to be run as administrator, but may yet fail as the license path has to be updated. This is done by selecting the “change” option instead of “repair” after running the executable file (install_quarc.exe), then pasting the updated license path. As of writing this, the location is: **C:\Global\Software\11920 QUARC Essentials Fh Vorarlberg.qlc**

Permission Denied

May occur if running the Simulink in a folder without sufficient permissions. Simply change current working directory to one you have permissions for :)